

# *AVChoice*

Veterinary Systems



*Universal*

# Analog systems ready for DR

## THE VETERINARIAN'S CHOICE

Universal AVChoice Systems deliver comprehensive veterinary imaging solutions designed for efficiency and cost-effectiveness.

## BUILT STRONG AND READY

Their sturdy welded construction tables can support patients up to 400 lbs. and feature smooth surfaces with urine trap perimeters

## COMPACT AND CLEAN

The AVChoice table cabinet can be used to house system generator components which greatly reduces floor space needed for installation. Complete system footprint is just 39" x 56" (with stationary or 2-way float top table option)

## HANDS-FREE OPERATION

Conveniently placed exposure switch and table top release strips are foot activated to ensure safety for the veterinarian and animals being examined. Release strips are located at both table ends and table front to maximize operability from all sides of the table.

Table lock release functions are electrically powered and quiet to prevent animals from being startled during examinations.

## STANDARD AVCHOICE SYSTEM

### TABLE

- Stationary Table Top

### TUBE STAND

- Fixed SID, 40" to image receptor
- Articulating arm on side of column for mounting generator control console

### GENERATOR (single or three phase)

- 30kW, high frequency, single phase, 400mA, 125kVp, with anatomical programming console
- 32kW, high frequency, three phase, 500mA, 125kVp, with anatomical programming console

### X-RAY TUBE

- 140kHu, 125kVp, 1.0-2.0mm focal spots

## SYSTEM OPTIONS

AVChoice Systems offer a variety of different component options to meet the needs of any veterinary practice.

### TABLE

- 2-way float top
- 4-way float top

### TUBE STAND

- Variable SID, 40" vertical range with electric locks

### GENERATOR

- 40kW and 50kW generator options available upon request
- Pedestal stand or wall mount for control console

### X-RAY TUBE

- 200kHu, 125kVp, 0.6-1.5mm focal spots

OPTIONAL – VARIABLE SID



ANGULATING TUBE ARM  
with angulation dial and operator handle



MANUAL COLLIMATOR  
with retractable tape measure

GENERATOR CONTROL CONSOLE  
Anatomically Programmed



ARTICULATING COLUMN MOUNT  
for generator control console


17" x 17" MANUAL CASSETTE TRAY  
for CR cassettes and cassette-sized DR detectors




EXPOSURE  
FOOTSWITCH

ELECTRIC TABLE TOP RELEASE STRIPS


**Front View**

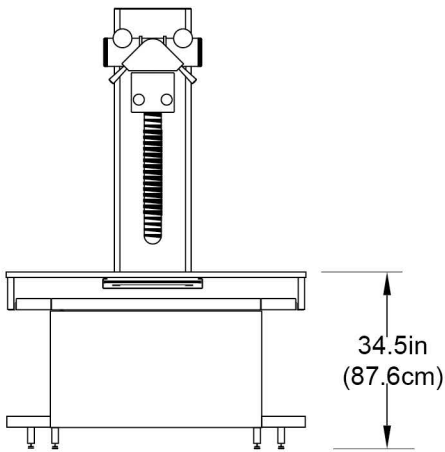


**Side View (Fixed SID)**

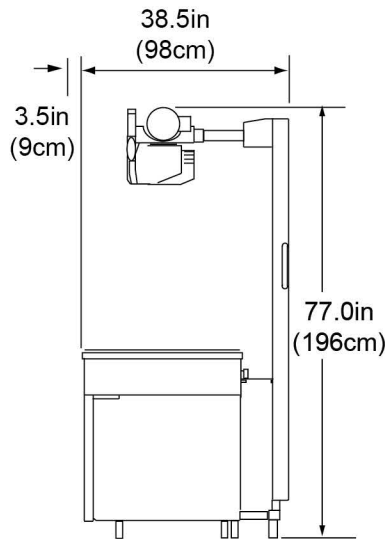

**Side View (Variable SID)**

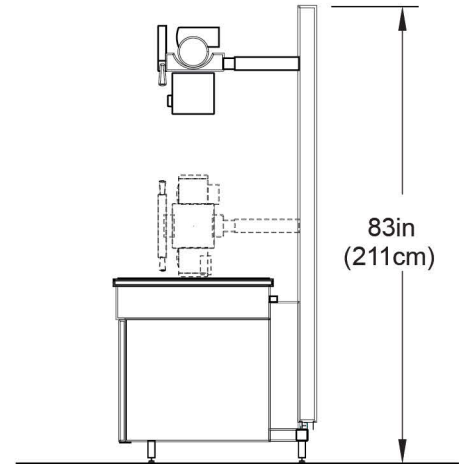

**Stationary and Two-way Float Top**

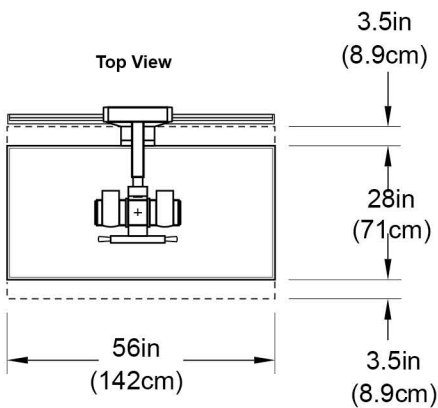

**Four-way Float Top**

